

Progranulin (Mouse) ELISA Kit

(Catalog # K4734-100; 100 assays; Store at 4°C)

I. **Introduction:** Progranulin (PGRN) also called epithelin precursor, proepithelin (PEPI), PC cell-derived growth factor (PCDGF), acrogranin, or paraganulin is a 593aa cysteine-rich protein of 68.5kDa, that is typically secreted in a highly glycosylated 88kDa form. BioVision's Progranulin (mouse) ELISA kit is a sandwich Enzyme Linked-Immunosorbent Assay (ELISA) for quantitative determination of mouse progranulin in biological fluids. A polyclonal antibody specific for progranulin has been precoated onto the 96-well microtiter plate. Standards and samples are pipetted into the wells for binding to the coated antibody. After extensive washing to remove unbound compounds, progranulin is recognized by the addition of a biotinylated polyclonal antibody specific for progranulin (Detection Antibody). After removal of excess biotinylated antibody, HRP labeled streptavidin (STREP-HRP) is added. Following a final washing, peroxidase activity is quantified using the substrate 3,3',5,5'-tetramethylbenzidine (TMB). The intensity of the color reaction is measured at 450 nm after acidification and is directly proportional to the concentration of progranulin in the samples. The lowest level of progranulin that can be detected by this assay is 60 pg/ml. The Assay range is 0.125 ng/ml – 8 ng/ml. This ELISA is specific for the measurement of natural and recombinant mouse progranulin. It does not cross-react with mouse RBP4, mouse adiponectin, mouse Nampt, mouse leptin, mouse RELM- β, mouse IL-33, mouse IDO, mouse resistin, mouse clusterin, mouse ANGPTL3, mouse ANGPTL4, human vaspin, human adiponectin, human progranulin, rat Nampt.

II. **Sample Type:** Serum and Cell culture supernatants.

III. **Kit Contents:**

Components	K4734-100	Part No.
Plate coated with mouse Progranulin antibody	6 x 16-wells strips	K4734-100-1
Wash Buffer (10x)	30 ml x 2	K4734-100-2
ELISA Buffer (10X)	30 ml x 2	K4734-100-3
Detection Antibody (DET)	20 μl	K4734-100-4
HRP Labeled Streptavidin (lyophilized)	2 μg	K4734-100-5
Mouse progranulin Standard (lyophilized)	16 ng	K4734-100-6
TMB Substrate Solution	12 ml	K4734-100-7
Stop Solution	12 ml	K4734-100-8
Plate sealers (plastic film)	2	K4734-100-9
Silica Gel Minibags	2	K4734-100-10


IV. **User Supplied Reagents and Equipment:**

- Microtiterplate reader at 450 nm, with the correction wavelength set at 540 nm or 570 nm
- Calibrated precision single and multi-channel pipettes. Disposable pipette tips
- Deionized water
- Microtubes or equivalent for preparing dilutions
- Disposable plastic containers for preparing working buffers
- Plate washer: automated or manual
- Glass or plastic tubes for diluting and aliquoting standard

V. **Storage and Handling:** Reagent must be stored at 2-8°C when not in use. Plate and reagents should be at room temperature before use. Do not expose reagents to temperatures greater than 25°C.

VI. **Reagent and Sample Preparation and Storage Conditions:**

- Bring all reagents and samples to room temperature (18 - 25°C) before use.
- Wash Buffer 10X has to be diluted with deionized water 1:10 before use (e.g. 50 ml Wash Buffer 10X + 450 ml water) to obtain Wash Buffer 1X.
- ELISA Buffer 10X has to be diluted with deionized water 1:10 before use (e.g. 20 ml ELISA Buffer 10X + 180 ml water) to obtain ELISA Buffer 1X.
- Detection Antibody (DET) has to be diluted to 1:1000 in ELISA Buffer 1X (10 μl DET + 10 ml ELISA Buffer 1X).
 NOTE: The diluted Detection Antibody is not stable and cannot be stored.
- HRP Labeled Streptavidin (STREP-HRP) has to be reconstituted with 100 μl of ELISA Buffer 1X.
 - After reconstitution of STREP-HRP, prepare aliquots and store them at -20°C. Avoid freeze/thaw cycles.
 - Dilute the reconstituted STREP-HRP to the working concentration by adding 8 μl in 12 ml of ELISA Buffer 1X (1:1500).
 - The diluted STREP-HRP is not stable and cannot be stored.
- Mouse Progranulin Standard (STD) has to be reconstituted with 1 ml of deionized water. This reconstitution produces a stock solution of 8 ng/ml. Mix the standard to ensure complete reconstitution and allow the standard to sit for a minimum of 15 minutes. Mix well prior to making dilutions. NOTE: The reconstituted standard is aliquoted and stored at -20°C
 - Dilute the standard protein concentrate (STD) (16 ng/ml) in Diluent 1X. A seven-point standard curve using 2-fold serial dilutions in Diluent 1X is recommended.
 - Suggested standard points are: 8, 4, 2, 1, 0.5, 0.25, 0.125, and 0 ng/ml.

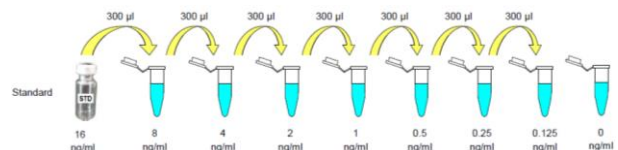

FOR RESEARCH USE ONLY! Not to be used on humans.

